


# Eye massager

## iQ100

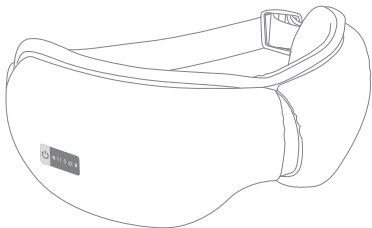

ambrint

---

[www.synergynext.co](http://www.synergynext.co)

# Product Highlights



5 Modes

3 Way massage system



Multi-zone massage

Double layered  
structure + 8 zone  
airbag massage



Thermostatic compress

Relieve eye fatigue



Built-in music

Bluetooth connects  
to a mobile phone  
to play music



One button operation

Easy control




Foldable body for  
on-the-go use

Easy to store Folding

# Contents of the box

Thank you for purchasing this product. Before use, please read the manual carefully to ensure that the machine, charging cable and accessories are complete and not damaged. Please pay special attention to the safety instructions.

Please fully charge the device before using it for the first time. It may take 3 to 4 hours.

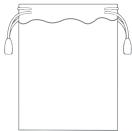

Ambrint Eye Massager Device

USB Type C Charging Cable


Travel Bag

User Manual

## How to wear

Before wearing the device on your eyes, please remove any eyewear, contact lenses, earbuds, hairpins or any other object which can create discomfort during the massage. After putting the device on, adjust the position of the device over your eyes and adjust the elastic band to a comfortable position.

Please wear the device over your eyes before turning it on, to prevent over-pressuring of the airbags.

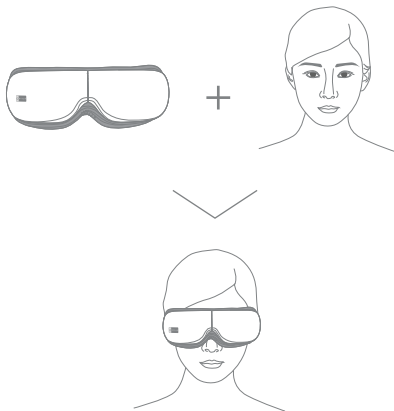

# How to use

## Switching on/off

- Press and hold the button for 3 seconds to switch on. The device will start in the default mode. A voice will inform you about the mode and the status of the device.
- Press and hold the button for 3 seconds to switch off.

## Mode switching

- Short press the button for 1 second to change the mode.
- Press again to select the next mode in sequence.

## Music

- Double click the button to turn off/on the built-in music.
- To listen to your own music, turn on the Bluetooth on your mobile phone (or any music device) and search for [YBY]. Connect to play. The device will start playing your own music.

## Auto off

- The device will automatically switch off after staying continuously on for 15 minutes.


## Charging

- Connect the USB Type C cable to the port and connect it to the charger.
- Red light is an indication of charging. It turns to blue when fully charged. During charging, the button function will be disabled and massage operatio will be unavailable. The device will not start if the charging cable is plugged-in to the charging port or the battery is too low.

# Selecting the mode

- Short press the button for 1 second to change the mode.
- Press again to select the next mode in sequence. A voice will inform you about the current mode.

(1. Integrated mode, 2. Clear mode, 3. Sleep mode, 4 Dynamic mode, 5. Comfort mode).

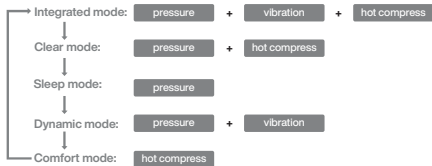

On/Off  
Mode selection  
Music control



The voice and the led lights will indicate the mode being used.



Bluetooth



Air pressure



Hot compress



vibration




Music pause/play

## Points to note

- This product is designed for personal use. It is not intended for medical or commercial purposes.
- People with special needs and disabilities should only use the device under special supervision.
- Children must be supervised so they do not use the device as a toy.
- It is recommended that you use the device for 15 minutes at a time.
- The device has a mild heat generating surface as a function of operation. People sensitive to heat must be cautious when using it.
- If you experience any physical abnormality during use, please stop using the device immediately and consult a doctor.
- Do not use the device in conjunction with other therapeutic devices.
- Do not use if you have eye problem, have had eye surgery, glaucoma, cataract, retinal detachment or other major eye diseases.
- Do not use if you have injuries around the eyes or on the forehead.
- Do not use the device if physically damaged by accident or by aging.

## Cleaning the device

- Clean the device surface with mild wet cotton towel to avoid scratches to the body and damage to the PU leather during cleaning.
- To prevent the device from makeup stains, and possible skin irritation, it is recommended to remove the face makeup before using the device.

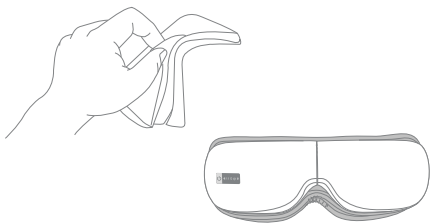

- Do not wash the device with water.
- Do not use cleaning fluids or chemicals to clean the device.



## Safety note

- Do not damage, twist, pull, knot or excessively bend the charging cable.
- Do not charge the device if it is likely to be submerged in water.
- Do not use a damaged power cable for charging. Do not attempt to repair damaged cable. Replace it with a standard Micro USB Type C cable.
- Store the device securely after daily use for product safety.
- To bring the device with you on Air travel, follow the airline instructions for placement in the luggage.
- It is recommended to keep the device in the hand luggage and not in the checked-in luggage for lithium battery air travel safety.
- Do not expose the device to excessive heat, fire, corrosive liquids or other adverse conditions.
- Save the manual and the carton box for future reference and storage.
- Store the device at a cool dry place.

# Product Specifications

## Ambrint Eye Massager

Model:	iQ100
Power:	3.4W
Battery:	Built-in
Battery:	Li-ion, 5V, 1200mAh