


Eye massager

iQ200

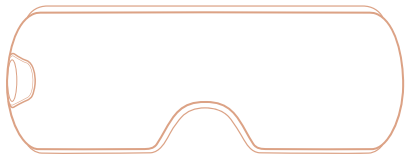

ambrint

Product Highlights



3D visualisation with
eye contour grooves



Multi-zone massage
Double layered
structure + 8 zone
airbag massage



Thermostatic compress
Relieve eye fatigue



Built-in music
Bluetooth connects
to a mobile phone
to play music




One button operation
Easy control



Foldable body for
on-the-go use
Easy to store Folding

Product construction

Thank you for using our products. Before use, please read the manual carefully to make sure that the machine, charging cable and other accessories are complete and undamaged. Please fully charge the device before using it for the first time.

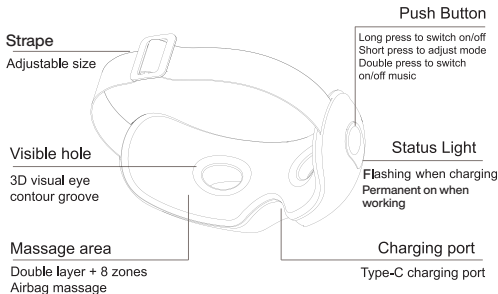

Type-C charging cable

How to use

Before wearing the device on your eyes, please remove any eyewear, contact lenses, earbuds, hairpins or any other object which can create discomfort during the massage. After putting the device on, adjust the elastic band to a comfortable position. Move the device up or down to align the see-through holes to the level of your eyes, enabling you to see without obstruction.

Please wear the device over your eyes before turning it on to prevent over-pressuring of the airbags.

Switching on/off

- Press and hold the button for 2 seconds to switch on. The device will start in the default mode. A female voice will inform you about the mode and the status of the device.
- Press and hold the button for 2 seconds to switch off.

Mode switching

- Short press the button for 1 second to change the mode.
- Press again to select the next mode in sequence. (Auto mode, Clear mode, Sleep mode).

Music

- Double click the button to turn off/on the built-in relaxing music.
- To listen to your own music, please turn on the Bluetooth on your mobile phone (or any music device) and search for [Eye Massager]. Connect to play, and the device will start playing your own music.

Smart Functions

- Low power indication

When the remaining power is left for about 10 minutes of operation, the status light will start flashing, along with an audible message, to remind you of recharging. The device will automatically switch off before the battery is fully exhausted.

- Charging indication

During charging, the status light will keep flashing, the button function will be disabled and massage operation will be unavailable. When the device is fully charged, the status light turns to always on.

- Automatic switch off

By default, the device automatically switches off after staying continuously on for 15 minutes.

Points to note

- Fully charge the device before using it for the first time. This may take up to 4 hours.
- This product is designed for personal use. It is not intended for medical or commercial purposes.
- People with special needs and disabilities should only use the device under special supervision.
- Children must be supervised to ensure they do not use the device as a toy.
- It is recommended that you use the device for 15 minutes at a time.
- The device has a mild heat generating surface as a function of operation. People sensitive to heat must be cautious when using it.
- If you experience any physical abnormality during use, please stop using the device immediately and consult a doctor.
- Do not use the device in conjunction with other therapeutic devices.
- Do not use if you have eye problem, have had eye surgery, glaucoma, cataract, retinal detachment or other major eye diseases.
- Do not use if you have injuries around the eyes or on the forehead.
- Do not use the device if it is physically damaged by accident or by aging.
- Do not charge the device if it is likely to be submerged in water.
- Do not damage, twist, pull, knot or excessively bend the power cable.
- Do not use a damaged power cable for charging. Do not attempt to repair a damaged cable. Replace it with a standard USB-C cable.

Safety Note

- Do not attempt to fully unfold the device, when the power cable is connected to the device charging port.
- Do not wash the device with water.
- Do not use cleaning fluids or chemicals to clean the device.
- To prevent harsh makeup stains, and possible skin irritation, it is recommended to remove the face makeup before using the device.
- Clean the device surface with mild wet cotton towel to avoid scratches to the body and damage to the PU leather during cleaning.
- Store the device in the pouch (provided with the product) after daily use for product safety.
- To bring the device with you on Air travel, follow the airline instructions for the luggage storage.
- It is recommended to keep the device in the hand luggage and not in the checked-in luggage for lithium battery air travel safety.
- Do not expose the device to excessive heat, fire, corrosive liquids or other adverse conditions.
- Save the carton box for future storage.
- Store the device at a cool dry place.

Product Specifications

Ambrint Eye Massager

Model:	iQ200
Power:	5W
Battery:	Built-in
Battery	Li-ion, 3.7V, 1200MAH

